

| | |
|---|---|
| EN | ECG ELECTRODES |
| DE | EKG-ELEKTRODEN |
| FR | ELECTRODES ECG |
| IT | ELETTRODI PER ECG |
| BG | ЕЛЕКТРОД ЗА ЕКГ |
| CS | EKG ELEKTRODY |
| DA | EKG-ELEKTRODER |
| EL | ΗΛΕΚΤΡΟΔΙΑ ΗΚΓ |
| ES | ELECTRODOS ECG |
| ET | EKG-ELEKTROODID |
| FI | EKG-ELEKTRODIT |
| HR | EKG ELEKTRODE |
| HU | EKG ELEKTRÓDÁK |
| KK | ЭКГ ЭЛЕКТРОДТАРЫ |
| LT | EKG ELEKTRODAI |
| LV | EKG ELEKTRODI |
| NL | ECG-ELEKTRODEN |
| PL | ELEKTRODY EKG |
| $PT_{EU}$ | ELÉTRODOS PARA ECG |
| $PT_{BRASIL}$ | ELETRODOS PARA ECG |
| RO | ELECTROZI EKG |
| SK | EKG ELEKTRÓDY |
| SL | EKG ELEKTRODE |
| SR | EKG ELEKTRODE |
| SV | EKG-ELEKTRODER |
| ID | ELEKTRODE EKG |
| JA | ECG電極 |
| KO | ECG 전극 |
| NO | EKG-ELEKTRODER |
| RU | ЭКГ ЭЛЕКТРОДЫ |
| TR | EKG ELEKTROTLARI |
| ZH | ECG 一次性使用心电电极 |
| UK | ЕЛЕКТРОДИ ДЛЯ ЕКГ |
| VI | ĐIỆN CỰC ĐIỆN TÂM ĐỒ |

GE HealthCare

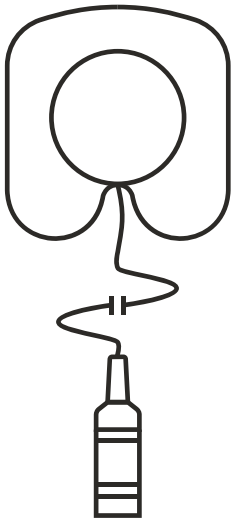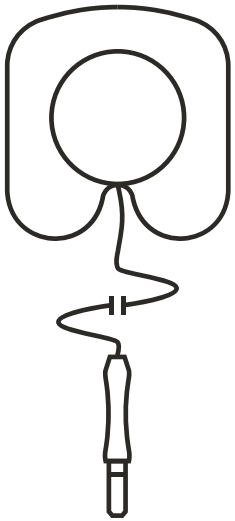Leonhard Lang GmbH, Archenweg 56, 6020 Innsbruck, Austria;
Fax +43 512 392210, medical@LeonhardLang.at

CE

93471 10076379_b 16.02.2023

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
## EN

**PRE-WIRED ECG ELECTRODES**
Please read this document and keep it safe. Ensure that any individual using this product is familiar with and understands the information contained in this document.

**INTENDED USE / INTENDED PURPOSE**
ECG electrodes are accessories for ECG monitoring devices. ECG electrodes are non-sterile, intended for single use on intact (uninjured) skin and are used to pick up ECG signals from the patient.

**WARNING**
- These pre-wired ECG electrodes must not be used during magnetic resonance imaging (MRI)! Non-compliance with this warning can lead to serious patient burns!
- To minimize the risk of burns at electrode placement sites during electrosurgery, always follow the procedure specified by the manufacturer of the electrosurgical or monitoring unit!
- Make sure that the patient does not become entangled in the cables! Beware of the risk of strangulation, especially with babies!

**CAUTION**
- Do not use once product has expired.
- Do not place the electrodes on areas of skin with inflammations, lesions, or injuries of any kind, or (for example) on skin stressed or damaged due to the use of medication!
- Do not use the electrodes with direct electrode contact on patients with a history of skin irritation, and report any incompatibility reactions to the specialist dealer.
- Do not prep the skin of patients with alcohol, solvents, or abrasive materials!
- Do not use electrodes if their gel has dried out!
- Do not use the electrodes longer than the time indicated on the packaging pouch!
- Only use each electrode once. In the event of reuse, the adhesive strength and electrical properties may be inadequate. There is also a risk of cross-infection from one patient to another.
- Use only electrodes that are suitable and appropriately labeled for X-ray imaging. These are marked with the “Radiolucent” symbol on the packaging pouch.

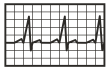**INSTRUCTIONS FOR USE**
The electrodes must only be used by a physician or a suitably trained member of medical staff.

**Selecting the Electrode Type:**
Use only electrodes that are suitable and appropriately labeled for X-ray imaging. These are marked with the “Radiolucent” symbol on the packaging pouch.

**Placement and Preparation:**
- Select the patient’s most favorable application site for the diagnostic or monitoring application in line with the hospital or medical specifications. Do not place electrodes on areas of the skin with inflammations, lesions, or injuries of any kind. Do not use the electrodes with direct electrode contact on patients with a history of skin irritation, and report any incompatibility reactions to the specialist dealer.
- Clean the selected application sites with water or with a soap suitable for the patient. Ensure that the application sites are dry and free of body oil in order to achieve optimum adhesion. Do not prep the skin of patients with alcohol, solvents, or abrasive materials.

**Applying the pre-wired ECG electrodes:**
- Remove the ECG electrodes from the pouch, disconnect the cables and straighten them.
- Remove the plastic film and apply the ECG electrodes to the selected positions. Gently smooth the electrodes with your finger as necessary in order to ensure full skin contact and optimum adhesion.
- Make sure that the patient does not become entangled in the cables! Beware of the risk of strangulation, especially with babies! The cables must not inhibit the patient’s extremities.
- Connect the electrode cables to the patient cable.
- Start the ECG.
- Always store any electrodes that are not immediately required in the unopened original pouch.

**Removing the pre-wired ECG electrodes:**
- After use, carefully remove the electrode using one hand whilst supporting the skin beneath it with the other. Loosen the electrodes at their corners or edges (not at the cable!) and slowly detach them. Ripping off the electrode or peeling it off quickly may damage the skin. Exercise particular caution if the skin is exceptionally sensitive, especially with children but also with elderly patients, diabetics, or patients on a prolonged course of certain medications known to cause dermatological side effects.
- Dispose of the product in accordance with local regulations or the guidelines of the hospital/health care facility.

For a patient/user/third party in the European Union and in countries with identical regulatory regime (Regulation 2017/745/EU on Medical Devices); if, during the use of this device or as a result of its use, a serious incident has occurred, please report it to the manufacturer and/or its authorized representative and to your national authority.

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